

Invacare® | Birdie® and Birdie® Compact



180kg weight limit
5 year warranty

A straightforward solution to safe and comfortable transfers

With maximum space for the client, both the *Invacare Birdie* and the *Birdie Compact* offer a truly comfortable lift and transfer to or from beds, chairs or even the floor. Ideal for use in domestic and residential care environments, the hoists are designed to fold and unfold easily without the need for tools.

Trusted Solutions, Passionate People





Comfort and ease of use combined

The high boom maximises space for the client even when the hoist is in the highest position. This gives an ergonomic lift and transfer as the client can be easily rotated 360°, which makes positioning much easier for the carer. The optimisation of space allows increased space in front of the actuator, reducing the risk of injury to the client's knees.

Ease of handling: simple to fold, unfold and dismantle without tools

Handling the hoist has been made as smooth and easy as possible as the hoists can easily be folded and unfolded without the need for tools. Hoists take up minimal space when stored and are easy to push and transport. If required, the mast and boom can be easily detached from the base (no tools required).

TÜV-certification for maximum security

For total security, the folding mechanism has a snap-lock system, thus minimizing the risk of accident. The **Birdie** and **Birdie Compact** are CE-certified in line with the Medical Device Directive 93/42/ECC.

To ensure the hoists comply with EN ISO 10535 (the European standard for hoists) the hoist has been rigorously tested both by Invacare and the German TÜV institute, an accredited external body.

Features and Options



Long boom - high user comfort

The client can easily be rotated and positioned. No risk of injuries to the knees.



Easy to fold and dismantle

Ease of transport and storage. No tools required when folding / unfolding or when dismantling the hoist.



Snap lock

Optimal security.



Emergency stop, lowering and integrated charger

Manual or electric emergency lowering (compact version only electric lowering). Control box is supplied with detachable or fixed battery.

Features and Options



Clip on spreader bar

Easy to change the spreader bar - just twist the safety catch, press to open - and then detach/attach the spreader bar.



Mounted spreader bar

Correctly and safe positioning of the spreader bar. A locking carabiner offers exceptional security.



Manual leg spread

Simply step on the pedal to adjust the leg spread.



Handle for leg spread

Handle can easily be attached to operate the manual leg spread.



Electrical leg spread

Available for *Invacare Birdie* models only.



2 point quick release spreader bars

Easy to attach and detach. Available widths: 35 cm, 45 cm, 55 cm.



Push bar and hand control

Superior ergonomics for carers.




Castors

100 mm castors (braking castor is shown).



Scale

If needed, a Homecare scale is available approved, according to medical standard EN 45501.

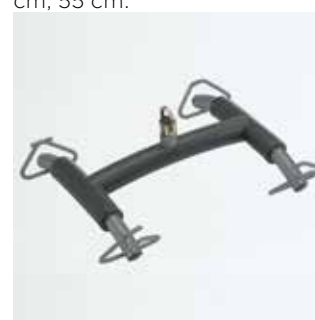

4 point quick release spreader bars

Easy to attach and detach. Available widths: 45 cm and 55 cm.

Technical data

| | | | | | |
|-----------------------|--------|-------|-------|-----------------|-------|
| | | | | | |
| Birdie | 124 cm | 64 cm | 42 kg | 44,5 - 171 cm | 10 cm |
| Birdie Compact | 109 cm | 52 cm | 36 kg | 52,5 - 161,5 cm | 10 cm |

| | | | | | |
|-----------------------|---------|--------|-----------|--------|--------|
| | | | | | |
| Birdie | 101 cm | 180 kg | 40 cycles | 2,9 Ah | 140 cm |
| Birdie Compact | 84,5 cm | 150 kg | 40 cycles | 2,9 Ah | 107 cm |

All measures are taken with 75 mm castors and 2 point spreaderbar. For 100 mm castors add 1,5 cm for height measures and 2 cm for width.

Frame colour



Light grey /
Dark grey

Warranty

5 year warranty:
Frame
Actuator
Control box

2 year warranty:
Hand control

Invacare Limited
Pencoed Technology Park -
Pencoed - Bridgend CF35 5AQ -
United Kingdom
Tel: +44 1656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: ireland@invacare.com
www.invacare.ie

© 2016 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Birdie - UK - 01/2016 - 1504431

The *Invacare Birdie* and *Birdie Compact* are part of the Safe Patient Handling range.

O'FLYNN MEDICAL

Authorised Distributer
O'Flynn Medical
The Dill Creamery, Tanyard,
Millstreet Town,
Co Cork, P51 D328.
T: 029 21789
Lo Call: 1890 440 440
E: info@oflynnmedical.com



www.oflynnmedical.com



Yes, you can.®