

"Kakadu" commode and shower chair made of aluminium



OY1. "Kakadu" commode and shower chair standard equipped with 4 pieces of 5" castors



OY2. "Kakadu" commode and shower chair equipped with 24" wheels with quick release (optional)

The "Kakadu" commode and shower chair has a frame made of white anodized aluminium and comes as standard with 4 pieces of 5" castors with parking brakes. The seat, backrest and armrests are made of soft black PU. The foot rests are easily adjustable in height and are removable.

All "Kakadu" chairs are supplied as standard with toilet bucket with carrying handle, lid and black PU seat cover.

The commode and shower chair can be equipped with 24" rear wheels with parking brakes. In the 24" wheel kit are included 2 brakes with clamps and brackets for fitting of the wheels on the frame. Pre-drilled holes makes the mounting of the wheels very easy. Anti-tip is optional.

Technical data	Small 302015	Standard 302018	Large 302019	XL 302020	Tilt 302024	Tilt L 302025	Diff. with 24"
Seat width	38 cm	45 cm	50 cm	55 cm	45 cm	50 cm	---
Seat depth	38 cm	43 cm	43 cm	43 cm	43 cm	43 cm	---
Seat height (adjust.)	50 - 55 cm	50 - 55 cm	50 - 55 cm	50 - 55 cm	55 - 60 cm	55 - 60 cm	---
Seat opening	17 x 33 cm	21,3 x 40,5	21,3 x 40,5	21,3 x 40,5	21,3 x 40,5	21,3 x 40,5	---
Pass over height	42 - 47 cm	42 - 47 cm	42 - 47 cm	40 - 45 cm	---	---	---
Pass over width	32 cm	40 cm	45 cm	50 cm	---	---	---
Total width **	49 cm	56 cm	61 cm	66 cm	56 cm	61 cm	+ 8 cm *
Total depth	85 cm	93 cm	93 cm	93 cm	110-120 cm	110-120 cm	+ 5 cm
Total height (adjust.)	96-101 cm	96-101 cm	96-101 cm	96-101 cm	97-141 cm	97-141 cm	---
Max. load	130 Kg	150 Kg	150 Kg	200 kg	130 Kg	130 Kg	---
Weight	11,00 Kg	12,90 Kg	13,50 Kg	14,10 kg	15,10 kg	15,70 kg	+ 2,90 kg
Tilt angle	---	---	---	---	-5° to +25°	-5° to +25°	---

* Total width: Difference for Kakadu Tilt (302024 or 302025) with 24" wheels mounted: + 14 cm

** +5 cm with the new accessory 30212M

"Kakadu tilt" commode and shower chair made of aluminium with tilt function



OY3. "Kakadu" commode and shower chair with tilt function



OY4. "Kakadu" commode and shower chair tilted

The commode and shower chair with tilt can be angled from -5° to + 25°. Head support with cut for ears, can be adjusted in height, depth and in the angle. The black padding is as the head support covered with soft PU. All "Kakadu" chairs are supplied as standard with toilet bucket with carrying handle, lid and black PU seat cover.

As standard the "Kakadu Tilt" is fitted with 5" wheels, but it can be equipped with 24" rear wheels with parking brakes. In the 24" wheel kit are included 2 brakes with clamps and brackets for fitting of the wheels on the frame. Pre-drilled holes makes the mounting of the wheels very easy. Anti-tip is optional.

Accessories for "Kakadu" commode and shower chair

For the commode and shower chair series "Kakadu" different accessories are available. See price list for the full range of accessories.



OY5. Amputation support




OY6. Security bar



OY7. Anti-tip
ATTENTION: usable only with mounted 24" wheels



OY8. Seat with round hole

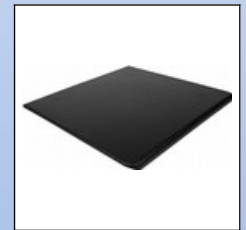

OY9. Seat cover for the whole seat

NEW "Commode Liner" - Hygiene bags for shower and toilet wheelchairs

The new patented **Commode Liner Disc** is a unique, portable, flexible shower chair accessory that is easy and hygienic to use.

The Liner Disc is used in conjunction with the **Commode Liner bags** (sold separately) and can be easily placed on and removed from the seat to prevent contamination of the seat and bucket.

The hygiene bags are available with plastic seat pad (picture OY10) or without (picture OY11).



OY10. Commode Liner Disc



OY11. Commode Liner