

SofCare® Chair Cushion Series

3-layer Technology



Convenient. Complete. Effective.

Stryker's SofCare family of products allows for a convenient, complete and effective way to enhance pressure redistribution and patient comfort while improving hospital outcomes and helping to prevent adverse events.

The SofCare Chair Cushion combines three layers of interlocking air cells to create a comfortably contoured support surface. SofCare elevates the body's sensitive pressure points off the surface of the chair and redistributes the patient's weight over the entire surface. The result is a superior combination of support and comfort.



The SofCare Chair Cushion has unique air-channeling technology which redistributes a patient's weight to reduce tissue interface pressure.



The SofCare Bariatric Chair Cushion has dual chambers for improved stability and optimal patient support. For nursing convenience, this cushion is also preinflated.


Pressure Management with 3-layer Technology

The unique 3-layer air-channeling technology creates interlocking air cells to immerse and envelop the patient providing low tissue interface pressure. The 3-layer design provides total air support preventing the patient from bottoming out because the patient is always supported by a pocket of air.



3-Layer Air-Channeling Technology

SofCare Chair Cushion Series

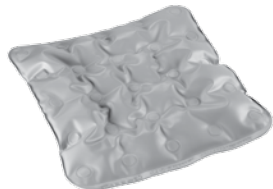

3800 E. Centre Ave.
Portage, MI 49002 U.S.A.
t: 269 389 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com

SofCare Chair Cushion Features

- Unique air-channelling technology redistributes a patient's weight to reduce tissue interface pressure
- Protects sensitive bony prominences susceptible to pressure ulcers
- Preinflated, ready to use
- Easy to transport
- Lightweight, waterproof and easy to clean
- Low-profile design fits all chairs
- Economical

SofCare Bariatric Chair Cushion Features

- Dual chambers for improved stability and optimal patient support
- Preinflated for nursing convenience
- Weight capacity up to 700 lb.
- Air adjustment valves for customization if needed.
- Small hand inflation pump available for pressure adjustments
- Same great features as the SofCare Chair Cushion
- Designed to fit most bariatric chairs

Specifications

| | SofCare Chair Cushion | SofCare Bariatric Chair Cushion |
|-----------------|-----------------------|---------------------------------|
| Model Number | CC476CE | CC2822 |
| Inflated Length | 18" (47 cm) | 22" (56 cm) |
| Inflated Width | 19.625" (50 cm) | 28" (71 cm) |
| Weight | 11.7 oz. (.33 kg) | 2.25 lb. (1 kg) |
| Weight Capacity | 350 lb. (158.8 kg) | 700 lb. (317.5 kg) |
| Packaged | 6 per carton | 5 per carton |
| Warranty | 30 days | 30 days |

Ordering Information

| | |
|---------|---------------------------------|
| CC476CE | SofCare Chair Cushion |
| CC2822 | SofCare Bariatric Chair Cushion |
| 91455 | Small Hand Inflation Pump |

Single Patient Use Only

Ask about our other SofCare 3-layer air-channelling technology products. Stryker reserves the right to change specifications without notice.