


SofCare System

Constant Low Pressure
Reactive Support Surface
Static Technology



Pressure Management with 3-layer Technology

The 3-layer, air-channeling technology creates interlocking air cells to immerse and envelop the patient providing low tissue interface pressure.

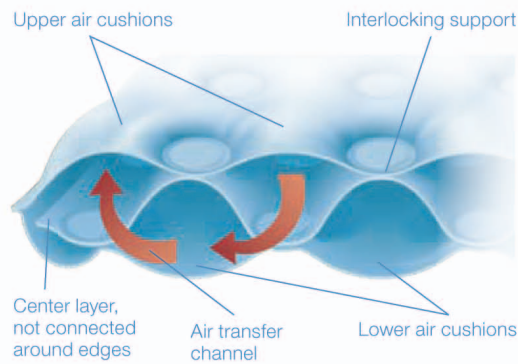

Constant Low Pressure

- Immersion
- Optimisation of the contact points
- No minimum weight
- Maximum internal pressure of 32 mmHg

3-layer Technology

300 Cells

Durable Monoblock Design



Control Unit

- Adjustable comfort/pressure
- Free flow air adjustment
- Max inflate
- Alarm (visible)

SofCare System



Cover Standard Features

- Four way stretched
- Breathable
- Waterproof

SofCare System

Model Number 2890 series

SofCare Specifications

Flame Retardant Standards	EN 597-1 and EN 597-2	
Safe Working Load	159 kg	(350 lb.)
Overall Length	188 cm	(74 in.)
Overall Width	89 cm	(35 in.)
Overall Height	15 cm	(5.9 in.)
Weight	5.00 kg	(11.0 lb.)

Pump Specifications

Power Supply	AC 230V 50 Hz, 0.05 A (for 230V system)	
Fuse Rating	T1AL, 250V	
Overall Length	25 cm	(9.8 in.)
Overall Width	12.5 cm	(4.9 in.)
Overall Height	8.5 cm	(3.3 in.)
Weight	1.22 kg	(2.7 lb.)

Classification

- Class IIa, Type BF, IP21
- Applied Part: Aire Mattress
- Not suitable for use in the presence of a flammable anesthetic mixture (No AP or APC protection)

Warranty

Two Years

*Stryker reserves the right to change specifications without notice.
Support surfaces are manufactured for Stryker Medical.*



O'Flynn Medical,
Millstreet, Co. Cork, Ireland.

Tel: 00353 (0)29 21799
Fax: 00353 (0)29 70191
Low Call: 1890 440 440
Email: info@oflynnmedical.com